



GTPB

The Gulbenkian Training Programme in Bioinformatics
(Since 1999)

Pedro Fernandes, Organiser



IB16S

Introductory Bioinformatics

12-16 December 2016

(Second 2016 run of this Course)

Basic Bioinformatics Sessions

Practical 4: Primer Design

Friday 9 December 2016

Primer Design

The prime intention of this exercise is to design a way to amplify a DNA fragment of reasonable size that includes a specific portion of the **PAX6** gene. The target region is that which includes the mutation you looked at earlier, that is a major cause of Aniridia. That is, the substitution that mutates the **33rd** amino acid position of the **PAX6** human protein. I remind that the details you discovered earlier are as follows:

<u>Affected Patient:</u>	<u>33rd amino acid position</u>	<u>mRNA Base position</u>
	Proline (P)	459
<u>Wild Type:</u>	<u>33rd amino acid position</u>	<u>Genomic DNA Base position</u>
	Alanine (A)	15714

The isolation, amplification and analysis of the target region of the genome could be affected by using **restriction enzymes**. In this case, there is more than one restriction enzyme whose cut site is dependant upon the mutation and so would produce a differing set of **restriction fragments** when used with the DNA of **Aniridia** affected patients to that normally expected. As long as those differences were course enough to be detected by a **Restriction Fragment Length Polymorphism (RFLP)** experiment. Software exists to select enzymes to isolate a chosen region of genomic DNA and to fragment that isolate in such a way it is possible to determine whether it includes the unfortunate mutation or not from the pattern of fragments generated.

For a variety of reasons, including the ready availability and ever decreasing cost of sequencing, this is typically not the preferred way to proceed. It is normally preferable to use **Polymerase Chain Reaction (PCR)** to isolate the region around the mutation and then to sequence samples from all individuals under examination. To do this, the first step would be to design suitable **PCR primers**. One program, in many different forms, is almost exclusively used for this purpose. The program is **primer3**. It is free and can be downloaded and run under linux and windows (at least). It is available as part of the **EMBOSS** package (**eprimer3**) and from a number of websites, including at the **Massachusetts Institute of Technology (MIT)**¹:

<http://bioinfo.ut.ee/primer3/>

This site is popular with many users wanting the very latest version of the software, complete control over the various options offered by **primer3** and are not too concerned with using a database search to check the uniqueness of the products they will produce.

Another excellent **primer3** web interface developed in the Netherlands is available at:


<http://www.bioinformatics.nl/cgi-bin/primer3plus/primer3plus.cgi>

The site incorporates access to a **blast** search to check the uniqueness of the selected primers (important if unwanted **PCR** products are to be avoided).

Mostly because of its completely seamless inclusion of a **blast** search to compare potential primers with appropriate sequence collections, I suggest we here use **primer3** as implemented at the **NCBI**, even though it offers less than complete control over the execution of **primer3** itself. **Go to:**

<http://www.ncbi.nlm.nih.gov>

Click on the **BLAST** option. Select

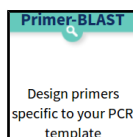

from the **Specialized BLAST** section.

¹ The MIT link here is to the latest version of **primer3** (version 4.0.0, soon **primer4** maybe?).

Upload your genomic **PAX6** sequence using the **Browse** (or **Choose File**) button for the **PCR Template**.

You have established that the mutation of greatest interest is the **G/C** substitution at position **15714** of the genomic sequence copied from **Ensembl**. It is logical therefore to specify that this feature be included in the **PCR** product not too near either end. Accordingly, request the **Forward primer** to be chosen **From** the region starting at base pair **15000** and continuing **To** base pair **15700**. Set the range for the **Reverse primer** to be **From 15800** and **To 16500**.

The default **PCR product size** is specified in the **Primer Parameters** section as between **70** and **1000** base pairs. This seems fine.

I would not presume to advise you on the melting temperatures that were most suitable². For this exercise, the defaults work splendidly.

By default, **primer-BLAST** will report the best **10** primer pairs it can find (**# of primers to return**). This is plenty for the exercise.

Do you think **10** primer pair suggestions is sufficient? If not, what number would you choose? _____

In addition to running **primer3** to suggest primers, **Primer-BLAST** checks against the possibility of unwanted **PCR** products by comparing potential primers against an appropriate sequence database with **blast**.

In the **Primer Pair Specificity Checking Parameters** section, set the **Database** selection to **RefSeq representative genomes**. Leave the **Organism** set as **Homo sapiens**.

You thus request each potential pair of **PCR** primers to be compared to the entire human genome. Thus unintended products of similar size to the intended product, can be identified.

The ideal conclusion is “just one product will be produced, on chromosome **11**, in the region of the **PAX6** gene”.


Use the appropriate ⓘ button to discover the purpose of the **Max target size** parameter.

This is a relatively new parameter replacing a very different parameter, the purpose of which was somewhat less obvious.

The reason for the **Max target size** parameter is surely pretty transparent, so maybe there is now less requirement to wake up its ⓘ button? For the present, the maximum size of any proposed **PCR** product, in this instance, is **1,000** base pairs (the form default). So the greatest size of an unwanted product that might be a problem (the **Max target size**) must be small enough to potentially be mistaken for a real product of **1,000** base pairs. **4,000** base pairs seems a bit cautious to me? However, unless you feel strongly about the matter, accept the default value of **4000**.

What value would you choose here if you were looking for uncluttered results? _____








2 My policy has been to not discuss parameters that pertain to the experimental conditions. In future versions of these notes, I will include discussion of some of these parameters. In the mean time, the ⓘ buttons are very helpful. I would also suggest the **MIT** site (or the **Wageningen** site) for very readable explanations linked from every parameter. The full **primer3** manual can be found [here](#).

Before setting **primer-BLAST** going, click on the **Advanced parameters** button. Not really so **Advanced**? More **Avoidable** by those in a hurry. At the top are the **Primer Pair Specificity Checking Parameters** that control the way that **blast** is run. Note the  buttons offering explanation.

Note the very high default **Blast expect (E) value**, suggesting you will be interested in matches with your primers that might occur up to **30000** times by chance! This does make sense as the primers will be very short and so many good, even exact, “chance” matches might be expected against a large database.

Comment upon the small default value for the **blast word size**? _____

Primer Pair Specificity Checking Parameters

Max number of Blast target sequences	50000	
Blast expect (E) value	30000	
Blast word size	7	
Max primer pairs to screen	500	
Max targets to show (for designing new primers)	20	
Max targets to show (for pre-designed primers)	1000	
Max targets per sequence	100	

Internal hybridization oligo parameters


☐ Pick internal hybridization oligo

	Min	Opt	Max
Hyb Oligo Size	18	20	27
Hyb Oligo tm	57.0	60.0	63.0
Hyb Oligo GC%	20.0	50	80.0

Note that you could get **primer-BLAST** to suggest an **Internal hybridisation oligo**, but decline the invitation this time.

Accept all the **Advanced parameters** as they are. Ask **primer-BLAST** to **Show results in a new window**.

Click on the **Get Primers** button.

Get Primers ☒ Show results in a new window ☒ Use new graphic view 

After a few moments of deep thought, **primer-BLAST** will notice that the template sequence you are using is **highly similar** (identical in fact) to part of an entry in the database being searched. Hardly surprising if one was to think about it.

You are invited to select all listed regions (just one this time) where matches with primers are likely to be the intended product. In this case, that is the whole list of one, so click on the **All** button. Every pair of primers that **primer3** selects **must** match this region

Input PCR template pax6-genomic sequence
Range 15000 - 16500

Your PCR template is highly similar to the following sequence(s) from the search database. To increase the chance of finding specific primers, please review the list below and select all sequences (within the given sequence ranges) that are intended or allowed targets.

Select: [All](#) [None](#) Selected: 1

Accession	Title	Identity	Alignment length	Seq. start	Seq. stop	Gene
<input checked="" type="checkbox"/> NC_000011.10	Homo sapiens chromosome 11, GRCh38.p2 Primary Assembly	100%	1501	31801962	31803462	PAX6


Submit ☒ Show results in a new window

of **Chromosome 11** as it is precisely the region investigated by **primer3** in the first place. This process avoids **blast** reporting intended products as unintended products.

Finally, all is ready, so ask to **Show results in a new window** and then click on the **Submit** button.

Once you have revelled in the opportunity to twiddle the fingers and scratch the ear(s) whilst **primers3** and **blast** go merrily about their appointed tasks, you will receive your results. These should look disarmingly like mine if all has gone well.

The summary **Graphic view** suggest just **3** solutions met the default criteria for success used by **primer3**. Up to **10** were permitted³.

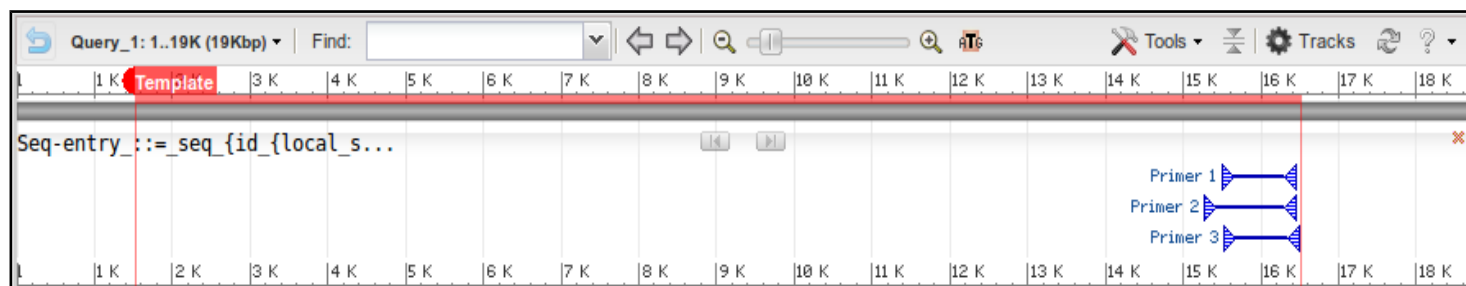

Hover your mouse over one or more and further details will pop up in separate windows.

Primer 1	Primer 2	Primer 3
<p>Forward: 15626..15645 length 20 Tm 59.75 GC 55.00% Seq AGGTCACAGCGGAGTGAATC</p> <p>Reverse: 16312..16331 length 20 Tm 58.77 GC 50.00% Seq TGACCTTGCTTAAAGTGCGC</p> <p>PCR product length: 706</p>	<p>Forward: 15409..15428 length 20 Tm 60.53 GC 60.00% Seq GATAGCAGGGAAGTACCGC</p> <p>Reverse: 16310..16329 length 20 Tm 58.31 GC 45.00% Seq ACCTTGCTTAAAGTGCGTT</p> <p>PCR product length: 921</p>	<p>Forward: 15651..15670 length 20 Tm 59.34 GC 55.00% Seq GGTGGTGTCTTTGTCAACGG</p> <p>Reverse: 16355..16376 length 22 Tm 59.84 GC 45.45% Seq GCCTGAAGATGCATCAAAACGA</p> <p>PCR product length: 726</p>

None of your three suggested primer pairs are reported with any unintended products, even given the very generous suggestion that products **4000** bases long should be considered a potential problem⁴.


Primer pair 1									
	Sequence (5'→3')	Template strand	Length	Start	Stop	Tm	GC%	Self complementarity	Self 3' complementarity
Forward primer	AGGTCACAGCGGAGTGAATC	Plus	20	15626	15645	59.75	55.00	6.00	3.00
Reverse primer	TGACCTTGCTTAAAGTGCGC	Minus	20	16331	16312	58.77	50.00	5.00	2.00
Product length 706									
Products on intended target									
>NC_000011.10 Homo sapiens chromosome 11, GRCh38.p7 Primary Assembly									
product length = 706									
Features associated with this product:									
paired box protein Pax-6 isoform a									
paired box protein Pax-6 isoform a									
Forward primer	1 AGGTCACAGCGGAGTGAATC	20							
Template	31802836	31802817							
Reverse primer	1 TGACCTTGCTTAAAGTGCGC	20							
Template	31802131	31802150							
Primer pair 2									
	Sequence (5'→3')	Template strand	Length	Start	Stop	Tm	GC%	Self complementarity	Self 3' complementarity
Forward primer	GATAGCAGGGAAGTACCGC	Plus	20	15409	15428	60.53	60.00	3.00	2.00
Reverse primer	ACCTTGCTTAAAGTGCGTT	Minus	20	16329	16310	58.31	45.00	5.00	1.00
Product length 921									
Products on intended target									
>NC_000011.10 Homo sapiens chromosome 11, GRCh38.p7 Primary Assembly									
product length = 921									
Features associated with this product:									
paired box protein Pax-6 isoform a									
paired box protein Pax-6 isoform a									
Forward primer	1 GATAGCAGGGAAGTACCGC	20							
Template	31803053	31803034							
Reverse primer	1 ACCTTGCTTAAAGTGCGTT	20							
Template	31802133	31802152							
Primer pair 3									
	Sequence (5'→3')	Template strand	Length	Start	Stop	Tm	GC%	Self complementarity	Self 3' complementarity
Forward primer	GGTGGTGTCTTTGTCAACGG	Plus	20	15651	15670	59.34	55.00	5.00	3.00
Reverse primer	GCCTGAAGATGCATCAAAACGA	Minus	22	16376	16355	59.84	45.45	8.00	2.00
Product length 726									
Products on intended target									
>NC_000011.10 Homo sapiens chromosome 11, GRCh38.p7 Primary Assembly									
product length = 726									
Features associated with this product:									
paired box protein Pax-6 isoform a									
paired box protein Pax-6 isoform a									
Forward primer	1 GGTGGTGTCTTTGTCAACGG	20							
Template	31802811	31802792							
Reverse primer	1 GCCTGAAGATGCATCAAAACGA	22							
Template	31802086	31802187							

- 3 Which rather makes mock of all the deep thought employed deciding upon the most sensible maximum number of predictions to be reported.
- 4 This was not true until very recently. **Primer-BLAST** reported many more primer pair suggestions and quite a few unintended products for each. The previous parameter restriction the length of unintended products was somewhat more generous.

As well as suggesting primers for PCR (or other purposes) and (optionally) suggesting hybridisation oligos, **primer-BLAST** can be used to evaluate user-selected primers. Earlier, you saved a pair of primer sequences associated with **PAX6** when searching the nucleotide databases at the **NCBI**. It would be interesting to discover the product these might produce. To do this you need an unsullied **Primer-BLAST** page. Go again to:

<http://www.ncbi.nlm.nih.gov>

Click on the **BLAST** option. Select

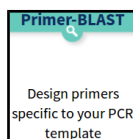

from the **Specialized BLAST** section.

Upload your genomic **PAX6** sequence using the **Browse** (or **Choose File**) button for the **PCR Template**.

Primer Parameters

Use my own forward primer (5'→3' on plus strand): CCAGCCAGAGCCAGCATGCAGAACA [Clear](#)

Use my own reverse primer (5'→3' on minus strand): GGTTGGTAGACTGGTGCTGAAACT [Clear](#)

PCR product size: Min 70 Max 1000

of primers to return: 10

Primer melting temperatures (T_m): Min 57.0 Opt 60.0 Max 63.0 Max T_m difference 3

Open up the file you made containing the primers from **GenBank** (**pax6_primers.fasta**) in a text editor.

Copy and Paste the two primer sequences into the **Use my own forward primer** and **Use my own reverse primer** boxes as appropriate.

Primer Pair Specificity Checking Parameters

Specificity check: ☒ Enable search for primer pairs specific to the intended PCR template

Search mode: Automatic

Database: Refseq representative genomes

Exclusion: ☐ Exclude predicted Refseq transcripts (accession with XM, XR prefix) ☐ Exclude uncultured/environmental sample sequences

Organism: Homo sapiens

Entrez query (optional):

Primer specificity stringency: Primer must have at least 2 total mismatches to unintended targets, including at least 2 mismatches within the last 5 bps at the 3' end. Ignore targets that have 6 or more mismatches to the primer.

Max target size: 20000

Splice variant handling: ☐ Allow primer to amplify mRNA splice variants (requires refseq mRNA sequence as PCR template input)

In the **Primer Pair Specificity Checking Parameters** section, set the **Database** selection to **RefSeq representative genomes**.


Leave the **Organism** as **Homo sapiens**.

Raise the **Max target size** parameter from **4000** to **20000**. You should check for enormous unintended products with this run of **Primer-BLAST**. The reasons for this will soon become apparent.

Ask primer-BLAST to Show results in a new window. Click on the **Get Primers** button.

Get Primers

☒ Show results in a new window ☒ Use new graphic view

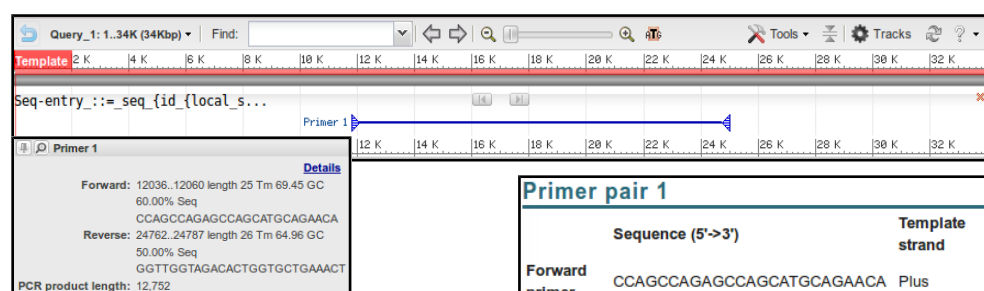

After a short thrill filled pause, you will receive a result that should again looks more that a trifle like mine.

Seemingly a fine match. Even the single **potentially unintended product** reported is actually the **intended product**. For some reason, **Primer-BLAST** does not eliminate intended products when searching for unintended ones when investigating user specified primers?

Success! However, applying a small measure of sober reflection, one has to wonder at a **PCR product** of **12,752** base pairs? I suspect that to be just a tad on the boastful side of probable? Clearly, **primer-BLAST** is convinced, but maybe a look at the references that came with these primer sequences would be advised before accepting this result at face value.

I have asked the guys at **NCBI** to explain. No full answer as yet, further prodding required. Prodded last **2016.04.02**. Maybe I give up?

Apparently, such a PCR product is possible! However, above **5,000** base pairs would be slow, require very close attention and be prone to errors.

Unfortunately, the only paper referenced does not explain what might be going on particularly clearly. However, there is a hint that the primers you saved were designed for use with mRNA/cDNA data. Therefore it might be interesting to run **primer-BLAST** one last time with **pax6_mrna.fasta** as the **PCR Template**.

Move back to your last **primer-BLAST** launch page. This time, load **pax6_mrna.fasta** as the **PCR Template**.

In the **Primer Pair Specificity Checking Parameters** section, set the **Database** selection set to **Refseq mRNA** and leave the organism set to **Homo sapiens**.

Set the **Max target size** back to its default value of **4000**, you should expect much smaller mRNA products this time, so no need for extending this maximum beyond **4000**.


These selections suppose that the design of **PCR** product was for selection from a library of all human cDNAs.

Ask **primer-BLAST** to **Show results in a new window**. Click on the **Get Primers** button.

Get Primers

Show results in a new window Use new graphic view

The result is a much more reasonable **Product length** of just **908** base pairs, reinforcing the theory that these primers were indeed designed for use with a cDNA library.

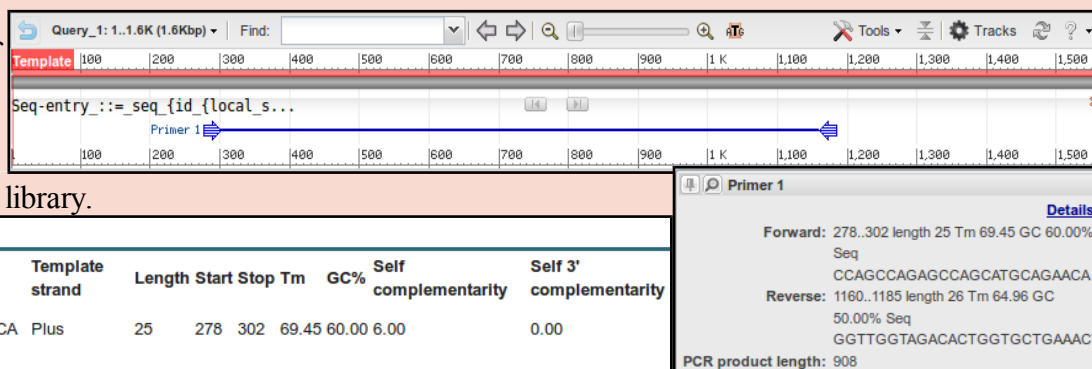

Primer pair 1									
	Sequence (5'→3')	Template strand	Length	Start	Stop	Tm	GC%	Self complementarity	Self 3' complementarity
Forward primer	CCAGCCAGAGCCAGCATGCAGAACA	Plus	25	278	302	69.45	60.00	6.00	0.00
Reverse primer	GGTTGGTAGACACTGGTGCTGAAACT	Minus	26	1185	1160	64.96	50.00	4.00	1.00
Product length	908								

Before moving on, afford a quick glance at the report offered concerning possible unintended products. Here **primer-BLAST** warns against human mRNAs that might be cloned along with the intended target.

Products on potentially unintended templates									
>NM_001310159.1 Homo sapiens paired box 6 (PAX6), transcript variant 9, mRNA									
product length = 908									
Forward primer	1	CCAGCCAGAGCCAGCATGCAGAACA	25						
Template	114	138						
Reverse primer	1	GGTTGGTAGACACTGGTGCTGAAACT	26						
Template	1021	996						
>NM_001310158.1 Homo sapiens paired box 6 (PAX6), transcript variant 8, mRNA									
product length = 950									
Forward primer	1	CCAGCCAGAGCCAGCATGCAGAACA	25						
Template	496	520						
Reverse primer	1	GGTTGGTAGACACTGGTGCTGAAACT	26						
Template	1445	1420						
>NM_001127612.1 Homo sapiens paired box 6 (PAX6), transcript variant 3, mRNA									
product length = 908									
Forward primer	1	CCAGCCAGAGCCAGCATGCAGAACA	25						
Template	455	479						
Reverse primer	1	GGTTGGTAGACACTGGTGCTGAAACT	26						
Template	1362	1337						
>NM_001258465.1 Homo sapiens paired box 6 (PAX6), transcript variant 7, mRNA									
product length = 908									
Forward primer	1	CCAGCCAGAGCCAGCATGCAGAACA	25						
Template	429	453						
Reverse primer	1	GGTTGGTAGACACTGGTGCTGAAACT	26						
Template	1336	1311						

The first thing to note is that the template (the mRNA sequence in the file **pax6_mrna.fasta**) is not a **RefSeq** mRNA. It comes from the **GenBank** database and so was included in the “non-redundant” union of databases you searched earlier.

Genbank sequences are generally generated directly from a specific sequencing project. **RefSeq** mRNAs are generally consensus sequences computed from the evidence represented by **Genbank** sequences. Consequently, there is no unintended product that we can ignore because it relates to the original template sequence.

All the unintended products could/would potentially be generated by the primers under investigation and have the potential to cause confusion. If you look down the list, you should conclude that the **9** unintended products come from **9** of the **11 RefSeq PAX6** transcripts found in the databases by test search and later detected by **blast**.

Why do you suppose **blast** did not pick up all the transcripts? _____

Note that the intended product is **908** base pairs long. Note that all the unintended products are either **908** long or **950** long. A difference of **42**.

How would you tell quickly which isoform was represented by each mRNA listed here? _____

Some fairly redundant questions to finish this section. I think I have already answered them all. But maybe you might wish to differ?

For all the “**potentially unintended products**”, the selected primers match exactly. Can you explain this? _____

```
>NM_001258463.1 Homo sapiens paired box 6 (PAX6), transcript variant 6, mRNA
product length = 908
Forward primer 1 CCAGCCAGAGCCAGCATGCAGAACA 25
Template 443 ..... 467
Reverse primer 1 GGTGGTAGACTGGTGTGCTGAACT 26
Template 1350 ..... 1325

>NM_001258463.1 Homo sapiens paired box 6 (PAX6), transcript variant 5, mRNA
product length = 950
Forward primer 1 CCAGCCAGAGCCAGCATGCAGAACA 25
Template 393 ..... 417
Reverse primer 1 GGTGGTAGACTGGTGTGCTGAACT 26
Template 1342 ..... 1317

>NM_001258462.1 Homo sapiens paired box 6 (PAX6), transcript variant 4, mRNA
product length = 950
Forward primer 1 CCAGCCAGAGCCAGCATGCAGAACA 25
Template 455 ..... 479
Reverse primer 1 GGTGGTAGACTGGTGTGCTGAACT 26
Template 1404 ..... 1379

>NM_001604.5 Homo sapiens paired box 6 (PAX6), transcript variant 2, mRNA
product length = 950
Forward primer 1 CCAGCCAGAGCCAGCATGCAGAACA 25
Template 443 ..... 467
Reverse primer 1 GGTGGTAGACTGGTGTGCTGAACT 26
Template 1392 ..... 1367

>NM_000280.4 Homo sapiens paired box 6 (PAX6), transcript variant 1, mRNA
product length = 908
Forward primer 1 CCAGCCAGAGCCAGCATGCAGAACA 25
Template 541 ..... 565
Reverse primer 1 GGTGGTAGACTGGTGTGCTGAACT 26
Template 1448 ..... 1423
```

DPJ – 2016.12.09

Model Answers to Questions in the Instructions Text.

Notes:

For the most part, these “**Model Answers**” just provide the reactions/solutions I hoped you would work out for yourselves. However, sometime I have tried to offer a bit more background and material for thought? Occasionally, I have rambled off into some rather self indulgent investigations that even I would not want to try and justify as pertinent to the objective of these exercises. I like to keep these meanders, as they help and entertain me, but I wish to warn you to only take regard of them if you are feeling particularly strong and have time to burn. Certainly not a good idea to indulge here during a time constrained course event!

Where things have got extreme, I am going to make two versions of the answer. One starting:

Summary:

Which has the answer with only a reasonably digestible volume of deep thought. Read this one.

The other will start:

Full Answer:

Beware of entering here! I do not hold back. Nothing complicated, but it will be long and full of pedantry.

This makes the Model answers section very big. **BUT**, it is not intended for printing or for reading serially, so I submit, being long and wordy does not matter. Feel free to disagree.

From your investigations of **Primer Design**

Do you think **10** primer pair suggestions is sufficient? If not, what number would you choose?

Until very recently, the default here was **5**. That seemed rather low to me. I included this question to solicit opinion rather than to impart knowledge. A default of **10** seems more in line with my instincts, but people who use this program seriously mostly tell me that they can select suitable primers from the first **2** or **3** suggestions of the program. So, **5** would seem a good choice and **10** would be moving towards cautiously overdoing things.

On the whole, informed opinion suggests that **10** suggestions will be more than enough in most circumstances.

What value would you choose here if you were looking for uncluttered results?

Summary:

Clearly, the smaller the number chosen, the shorter will be the list of spurious products. However, pick something too small and you risk including unintended product(s) that could cause confusion. The size selected must be sufficient that larger unwanted PCR product(s) could easily be spotted by other means (simply by size?).

Full Answer:

Well, mostly for me, and just in case you were curious, when I first wrote the question, the parameter was very different and not so easy to understand. Pure self indulgence, I know, but here is the history. The parameter explained itself, via the ⓘ button, thus:

Misprimed product size deviation	<input type="text" value="4000"/> ⓘ This specifies the size variation of the off-target PCR products relative to that of your intended PCR product. Only those primer pairs producing an off-target PCR product within the specified range will be tagged as non-specific.
----------------------------------	---

I interpreted this to mean that only **blast** predicted products of up to **X+4,000** base pairs, where **X** base pairs is the length of the intended target, will be given any regard. It is thus assumed that a difference of **4,000** base pairs between an intended PCR product (predicted by **primer3**) and a spurious product (detected by **blast**) can easily be detected simply by size difference.

Of course this parameter also will reject unwanted **blast** predicted products that are less than **X-4,000** base pairs will be given any regard. Given the largest possible **primer3** suggestion will be **1,000** base pairs (the form setting for the exercise specifies products of between **100**⁷ and **1,000** base pairs), this is hardly an issue here.

Comment upon the small default value for the **blast word size**?

By default, **blast** will be looking for aligned exactly matching blocks of **7** nucleotides when identifying where a primer might match a database entry. The entire primer match with the template sequence does not have to be exact for the primer to be acceptable. The entire primer is typically only around **20** bases long. And word size much more than **7** would clearly miss too much to be effective.

Why do you suppose **blast** did not pick up all the transcripts?

Well, the simple answer is that the transcripts that were not detected as unwanted products cannot include either the forward primer, or the reverse primer, or both. This is, almost, the only possible explanation.

I have a huge **Full Answer:** for this one that needs either throwing away or bringing up to date. I think I was the only one who ever read it, so you miss little while I decide its fate.

How would you tell quickly which isoform was represented by each mRNA listed here?

All the mRNAs reported were of length **908** or **950**.

A reasonable guess might be based on the length of the products? All those that are **908** bases might be assume to produce the **422** amino acid **canonical isoform**. All those that are **950** (i.e. **42** base pairs longer) might be assumed to **436** produce amino acid **isoform 5a** proteins (i.e. **14** amino acids longer).

Just a guess of course, but one I would be happy to have faith in. To be certain, one would need to read the annotations of each listed **RefSeq** entry!

⁷ The form explicitly declares a minimum of **70**, but the ranges from which the **forward** & **reverse** primers must come (**15000-15700** & **15800-16500**) make the smallest possible **primer3** prediction **100** base pairs long.

Is the number of “**potentially unintended products**” as you would expect, given the evidence from **Entrez** and **blast**?

Yes, I think so, given you accept my investigation (see above) as to why there were only **16** “**potentially unintended products**” when you might have expected **24**, given your **blast** results. **GeneCards** now encourages an initial expectation of **24** “**potentially unintended products**”. **Ensembl** only uses the higher quality **RefSeq** mRNAs. Currently, **Ensembl** uses **10** of the **11** good quality **RefSeq** mRNAs to make its transcripts predictions. Close enough?

For all the “**potentially unintended products**”, the selected primers match exactly. Can you explain this?

Well, of course they do??? All the transcripts found are generated from the same region of genomic DNA and therefore will be identical in all shared regions, including the primer regions. I suppose, in other instances, it would be possible to have transcripts with variation in the regions matching the primers insufficient to stop the primers working? But not in this case.

One might conclude there are no genuinely “unintended” products? All are real **PAX6** transcripts of varying certainty. A genuine unintended product would come from an entirely different part of the genome and would not necessarily match exactly with respect to the primers. They would just need to be “good enough to work”.

The “**potentially unintended products**” are of different sizes. Can you explain the difference between the possible product lengths?

Are the numbers of “**potentially unintended products**” of each possible length consistent with your **blast** results?

Yes yes yes! I think both these questions made a bit more sense a few generations of these notes ago. We have already answered them sufficiently I suggest. I refer you to the answers above.

DPJ – 2016.12.09